MPFD Foil Activation Experiment Resource

John Boyington Kansas State University

Inventory

- 1. Loaded wand
- 2. Stopwatch
- 3. Sample Bags (Four plastic bags labeled #1, #2, #3, and #4; #1 corresponds to the highest axial position.)
- 4. Sample pig for transportation
- 5. Wire Cutters

Procedures


For each wand:

- 1. Raise reactor power to P then insert wand into reactor core.
- 2. Irradiate foil for t_i .
- 3. Scram reactor and remove wand from reactor core.
- 4. Place in fuel storage rack for t_w .
- 5. Pull wand to surface and remove internals. Remove foils from internals and place in labeled bags. Place bags in sample pig.
- 6. Count foils one at a time with HPGe for listed counting time.

Irradiation Order (6/22/18)

- 1. Wand 0 Rh, Au, In
- 2. Wand 2 In(Cd)
- 3. Wand 3 Au(Cd)
- 4. Wand 4 Al

Principle Reactions

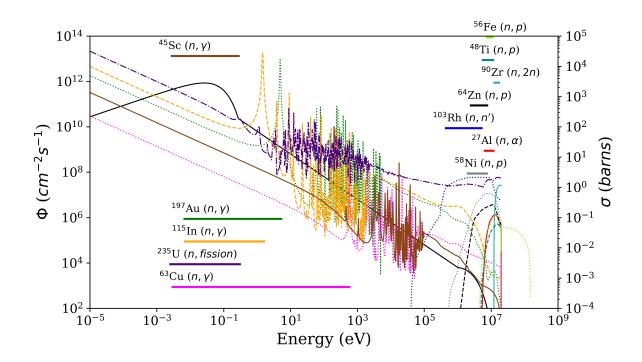

| Reaction | $T_{1/2}$ | ROI (eV) | Important Gammas (keV) |
|--|-----------|------------------------|----------------------------------|
| $^{115} In(n,\gamma)^{116} In$ | 54 m | 7.0021e-03, 1.6130e+00 | 417, 819, 1090, 1293, 1508, 2111 |
| 115 In(n, γ) 116 In Cd | 54 m | 1.1955e+00, 1.9916e+00 | 417, 819, 1090, 1293, 1508, 2111 |
| $^{197}\mathrm{Au}(\mathrm{n},\gamma)^{198}\mathrm{Au}$ | 2.7 d | 6.7266e-03, 5.2684e+00 | 412, 676, 1088 |
| $^{197}\mathrm{Au}(\mathrm{n},\gamma)^{198}\mathrm{Au}\;\mathrm{Cd}$ | 2.7 d | 4.0752e+00, 7.1730e+00 | 412, 676, 1088 |
| 103 Rh(n,n') 103m Rh | 56.12 m | 4.4469e+05, 5.1947e+06 | 40 |
| $^{27}\mathrm{Al}(\mathrm{n},\alpha)^{24}\mathrm{Na}$ | 15.03 h | 6.4564e+06, 1.1695e+07 | 1369, 2754 |


Indium

Power Level: 100.0 kW(th)

Time at Power: 60 s Wait Time: 3600 s

Total Activity at Removal: 1.79e+04 μCi

| Position | Mass mg | Start Counting s | Counting Time s | Counting Activity μCi |
|----------|-----------|------------------|-------------------|----------------------------|
| 1 | 1.7 | 3660 | 300 | 2.27e + 03 |
| 2 | 1.5 | 3960 | 300 | 1.88e+03 |
| 3 | 1.4 | 4260 | 300 | 1.65e + 03 |
| 4 | 1.6 | 4560 | 300 | 1.77e + 03 |

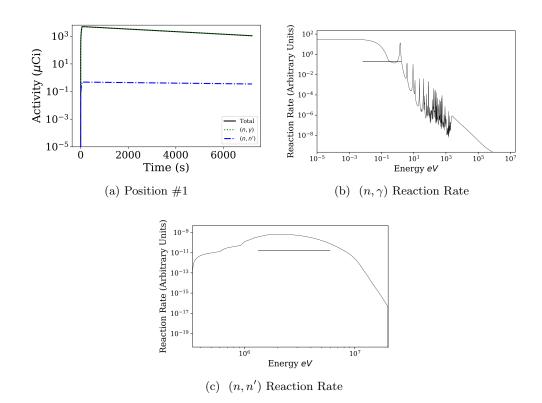

| Reaction | $T_{1/2}$ | ROI (eV) | Important Gammas (keV) |
|---------------|-----------|--------------------|---|
| (n, γ) | 54.0 m | 7.00e-03, 1.61e+00 | 138(0.03), 417(0.36), 819(0.17), 1090(0.53), 1293(0.8), 1508(0.11), |
| (n, n') | 4.4 h | 1.33e+06, 5.96e+06 | 335(0.5) |


Indium (Cd)

Power Level: 100.0 kW(th)

Time at Power: 60 s Wait Time: 3600 s

Total Activity at Removal: 8.83e+03 μCi

| Position | Mass mg | Start Counting s | Counting Time s | Counting Activity μCi |
|----------|-----------|------------------|-------------------|----------------------------|
| 1 | 1.7 | 3660 | 300 | 1.12e+03 |
| 2 | 1.5 | 3960 | 300 | 9.30e+02 |
| 3 | 1.4 | 4260 | 300 | 8.14e+02 |
| 4 | 1.6 | 4560 | 300 | 8.73e + 02 |

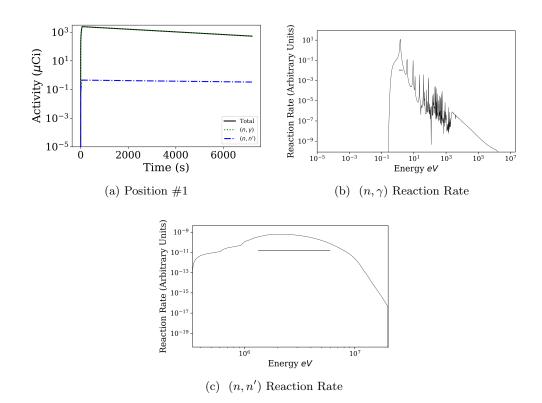


| Reaction | $T_{1/2}$ | ROI (eV) | Important Gammas (keV) |
|---------------|-----------|--------------------|---|
| (n, γ) | 54.0 m | 1.20e+00, 1.99e+00 | 138(0.03), 417(0.36), 819(0.17), 1090(0.53), 1293(0.8), 1508(0.11), |
| (n, n') | 4.4 h | 1.34e+06, 5.97e+06 | 335(0.5) |

Gold

Power Level: 100 kW(th) Time at Power: 60 s Wait Time: 3600 s

Total Activity at Removal: 2.78e+03 μCi

| Position | Mass mg | Start Counting s | Counting Time s | Counting Activity μCi |
|----------|-----------|------------------|-------------------|----------------------------|
| 1 | 5.0 | 3660 | 300 | 5.60e+01 |
| 2 | 4.35 | 3960 | 300 | 4.87e + 01 |
| 3 | 4.3 | 4260 | 300 | 4.81e+01 |
| 4 | 4.37 | 4560 | 300 | 4.88e + 01 |

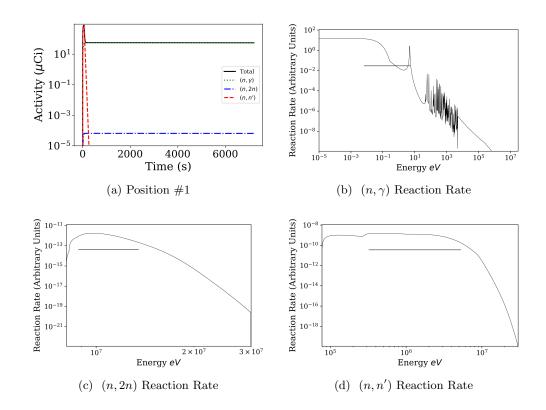


| Reaction | $T_{1/2}$ | ROI (eV) | Important Gammas (keV) |
|--------------|-----------|--------------------|--|
| (n,γ) | 2.7 d | 6.73e-03, 5.27e+00 | 412(0.95), 676(0.01), 1088(0.002) |
| (n,2n) | 6.2 d | 8.83e+06, 1.35e+07 | 333(0.25), 356(0.94), 426(0.06), 1091(0.002) |
| (n, n') | 7.8 s | 3.24e+05, 5.25e+06 | 130(0.08), 279(0.75) |

Gold (Cd)

Power Level: 100 kW(th) Time at Power: 60 s Wait Time: 3540 s

Total Activity at Removal: 2.52e+03 μCi

| Position | Mass mg | Start Counting s | Counting Time s | Counting Activity μCi |
|----------|-----------|------------------|-------------------|----------------------------|
| 1 | 5.0 | 3600 | 300 | 3.01e+01 |
| 2 | 4.35 | 3900 | 300 | 2.62e+01 |
| 3 | 4.3 | 4200 | 300 | 2.59e+01 |
| 4 | 4.37 | 4500 | 300 | 2.63e+01 |

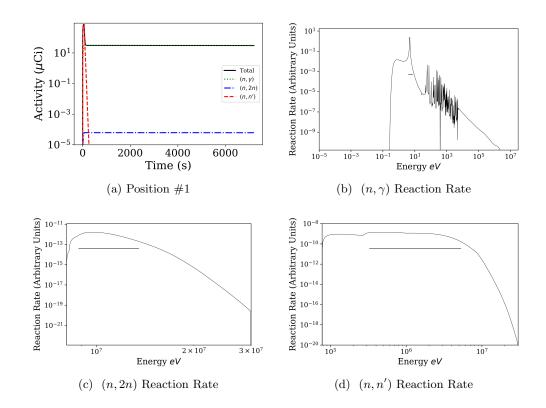


| Reaction | $T_{1/2}$ | ROI (eV) | Important Gammas (keV) |
|--------------|-----------|--------------------|--|
| (n,γ) | 2.7 d | 4.08e+00, 7.17e+00 | 412(0.95), 676(0.01), 1088(0.002) |
| (n,2n) | 6.2 d | 8.83e+06, 1.35e+07 | 333(0.25), 356(0.94), 426(0.06), 1091(0.002) |
| (n, n') | 7.8 s | 3.28e+05, 5.28e+06 | 130(0.08), 279(0.75) |

Rhodium

Power Level: 100 kW(th) Time at Power: 60 s Wait Time: 3600 s

Total Activity at Removal: 4.02e+04 μCi

| Position | Mass mg | Start Counting s | Counting Time s | Counting Activity μCi |
|----------|-----------|------------------|-------------------|----------------------------|
| 1 | 0.7 | 3660 | 600 | 1.95e+00 |
| 2 | 0.55 | 4260 | 600 | 8.41e-01 |
| 3 | 0.5 | 4860 | 600 | 5.79e-01 |
| 4 | 0.55 | 5460 | 600 | 5.41e-01 |

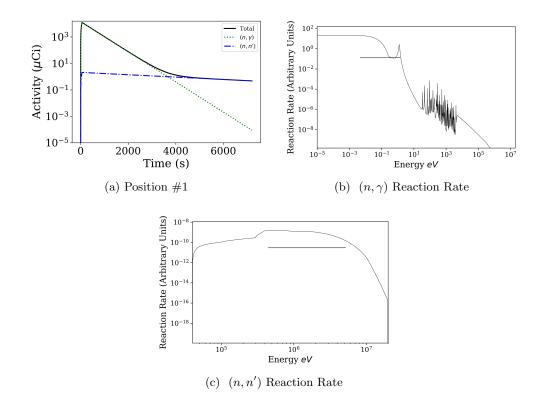


| Reaction | $T_{1/2}$ | ROI (eV) | Important Gammas (keV) |
|--------------|-----------|--------------------|--|
| (n,γ) | 4.4 m | 4.67e-03, 1.39e+00 | 51(0.47), 78(0.025), 560(0.026), 770(0.0018) |
| (n, n') | 56.1 m | 4.45e+05, 5.19e+06 | 40(0.004) |

Aluminum

Power Level: 100 kW(th) Time at Power: 3600 s Wait Time: 234000 s

Total Activity at Removal: $3.79e+02~\mu Ci$

| Position | Mass mg | Start Counting s | Counting Time s | Counting Activity μCi |
|----------|-----------|------------------|-------------------|----------------------------|
| 1 | 0.3 | 237600 | 3600 | 2.07e-04 |
| 2 | 0.2 | 241200 | 3600 | 1.32e-04 |
| 3 | 0.1 | 244800 | 3600 | 6.30e-05 |
| 4 | 0.2 | 248400 | 3600 | 1.20e-04 |

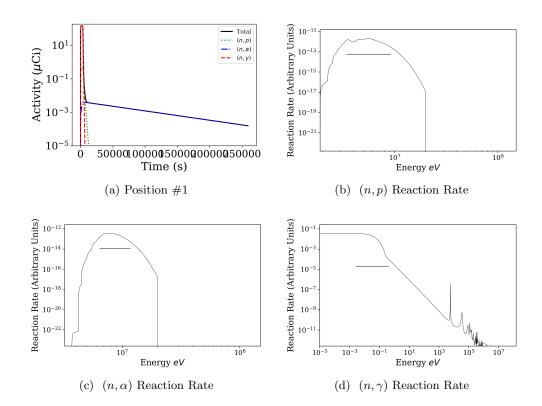

| Reaction | $T_{1/2}$ | ROI (eV) | Important Gammas (keV) |
|---------------|-----------|--------------------|---------------------------------|
| (n,p) | 9.5 m | 3.42e+06, 9.14e+06 | 180(0.007), 840(0.7), 1013(0.3) |
| (n, α) | 15.0 h | 6.46e+06, 1.17e+07 | 1369(1), 2754(1) |
| (n, γ) | 2.2 m | 2.82e-03, 4.15e-01 | 1780(1) |

Useful Links

Activation Calculator

https://www.ncnr.nist.gov/resources/activation/

Online Spectrum Catalogs for Ge and Si(Li)

http://www4vip.inl.gov/gammaray/catalogs/ge/catalog_ge.shtml

Decay Radiation Search

https://www.nndc.bnl.gov/nudat2/indx_dec.jsp

Evaluated Nuclear Data File (ENDF) Retrieval & Plotting https://www.nndc.bnl.gov/sigma/